

Product Name:	50mL Mini Bioreactor		
Product Application:	This product is a test tube specially used to grow cells. It has good chemical resistance and high mechanical strength. Suitable for any 50ml standard centrifugal device.		
Product Materials:	HDPE (cap) + PP (tube) Meets <i>USP, Class VI standards</i>		
Product Shelf Life:	Temperature range: stored at room temperature Shelf life: 3 years after date of production		
Sterilization:	Yes. Sterilized by E-beam, Sterility Assurance Level: SAL=10 ⁻⁶ . The product has been irradiated and dosimetrically released based on ANSI/AAMI/ISO 11137		
Pyrogens:	Non-Pyrogenic		
RNase/DNase:	DNase/RNase free		
BSE/TSE:	These products are deemed animal free.		
Performance Testing:	Each manufacturing lot is sampled and tested in accordance with standard operating procedures. Appearance inspection: qualified Sealing test: qualified Cell culture test: qualified Packaging inspection: qualified		
Features	1. High transparency, easy to observe. 2. Vent cap with 0.22 µm hydrophobic filter membrane which is anti-bacterial and breathable. 3. High volume dissolved oxygen coefficient and fast gas exchange. 4. Can be used for centrifugation, RCF 12000*g 5. Withstand temperature range: Approx.-80°C to 60°C. 6. Packaged in a sterile ziplock bag. 7. Product batch number identification for easy quality traceability		

Head Office

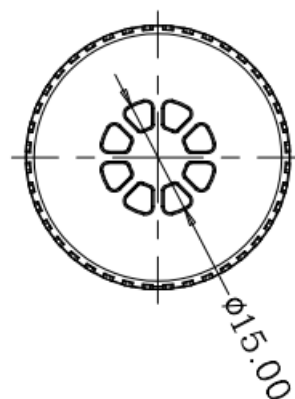
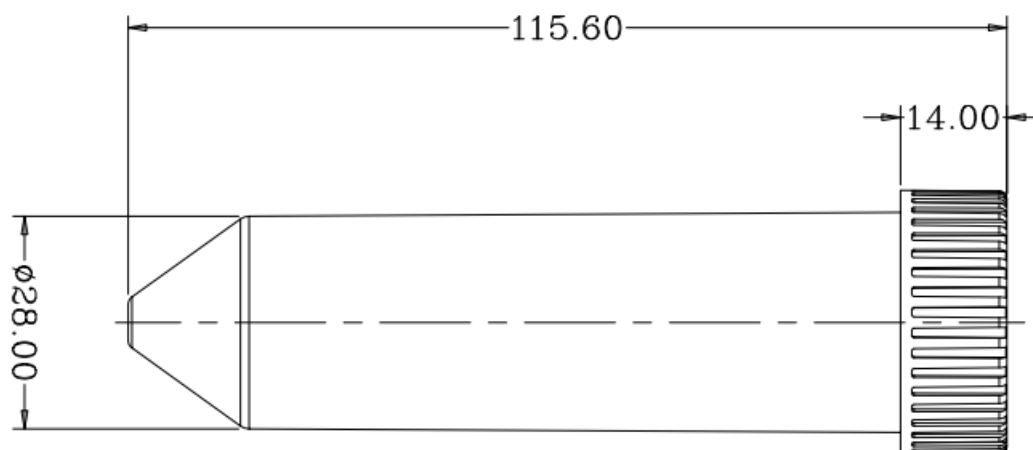
Email: info@cell-nest.com
Online: www.cell-nest.com

Overseas

NEST USA (New Jersey/ Phoenix)
NEST Scientific 株式会社 (Yokohama, Japan)
NEST Scientific Europe B.V (Netherlands)
NEST Scientific (MENA) FZE (Sharjah, United Arab Emirates)

Product Range

Cat. No.	Volume (mL)	Packaging	Max.RCF(xg)	/Pack	/Case
788211	50	Bagged	12000	10	100

Technical Drawing**Head Office**

Email: info@cell-nest.com
Online: www.cell-nest.com

Overseas

NEST USA (New Jersey/ Phoenix)
NEST Scientific 株式会社 (Yokohama, Japan)
NEST Scientific Europe B.V (Netherlands)
NEST Scientific (MENA) FZE (Sharjah, United Arab Emirates)